




## **SULFO CY3 NHS ESTER**

**SKU:** FP-1656

## **Description**

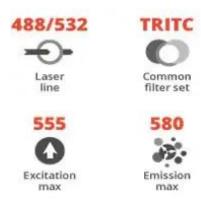

Sulfo-Cy3 NHS Ester is a highly hydrophilic bright, and pH insensitive from pH 4 to pH 10 orange-fluorescent dye. Sulfo-Cy3 conjugates of antibodies, peptides, and proteins can be excited using the 532 nm or 555 nm laser line and visualized with TRITC (tetramethylrhodamine) filter sets. Sulfo-Cy3 conjugates give less background than TAMRA and most other commonly used fluorescent dyes.

The NHS ester (or succinimidyl ester) is the most popular amine reactive group for labeling with the primary amines of proteins (Lys), amine-modified oligonucleotides, and other aminecontaining molecules. Sulfo-Cy3 NHS ester is not recommended for labeling proteins at high

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.

Telephone: (650) 697-3600








molar ration due significant self-quenching, and only recommended for detection of moderate-to-high abundance targets. For detection of low-abundance biological targets we recommend using AZDye 555 NHS Ester (Alexa Fluor ® 555 analog), which can be attached to proteins at high molar ratios without significant self-quenching, enabling brighter conjugates and more sensitive detection.

## Abs/Em Spectra

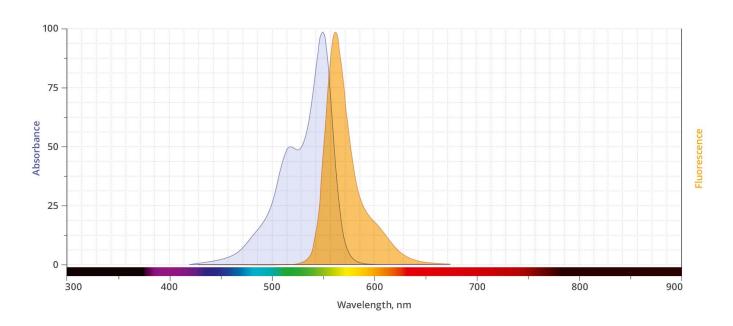

## **Specifications**

**Unit Size** 5 mg, 25 mg, 100 mg

**Reactivity** Primary amine

**Abs/Em Maxima** 555/572 nm

**Extinction coefficient** 150,000 cm-1M-1 **Solubility** Water, DMSO, DMF

Spectrally similar dyes Cy3, DyLight® 555, Alexa Fluor® 555

Molecular weight 821.93 (protonated)

**Storage Conditions** -20°C.

**Shipping Conditions** Ambient temperature

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.





Email: customerservice@vectorlabs.com

Telephone: (650) 697-3600

For research use only. Not intended for therapeutic or diagnostic use in animals or humans.