


## **VECTASHIELD® ANTIFADE MOUNTING MEDIUM**

**SKU:** H-1000-10

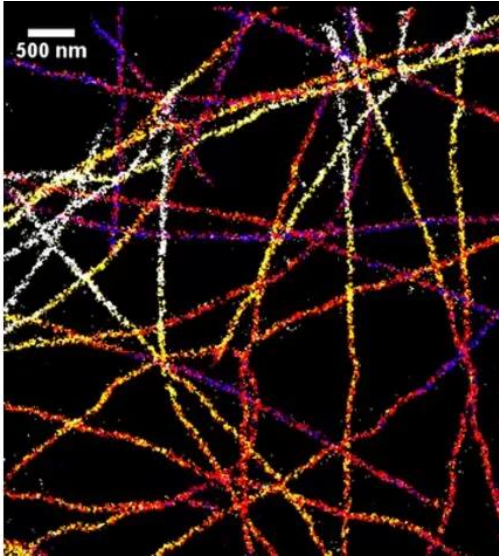

### **Description**

VECTASHIELD Antifade Mounting Medium is a unique, stable formula for preserving fluorescence. This medium inhibits rapid photobleaching of fluorescent proteins and fluorescent dyes.

### **Features:**

- Ideal refractive index
- Ready to use
- No warming necessary
- Can be stored without sealing for long term analysis
- Easy to use

The original VECTASHIELD Mounting Medium does not solidify, but remains a liquid on the slide and can be stored without sealing. If desired, coverslips can be sealed around the perimeter with nail polish or a plastic sealant. Mounted slides should be stored at 4 °C, protected from light.


**For research use only. Not intended for therapeutic or diagnostic use in animals or humans.**

## Specifications

<b>Unit Size</b>	10 ml
<b>Applications</b>	Immunofluorescence, In situ hybridization, Cellular Imaging, Super-Resolution Microscopy
<b>Mounting</b>	Aqueous (Non-Hardening)
<b>Antifade</b>	Yes
<b>Counterstain</b>	None

---

## GALLERY IMAGES

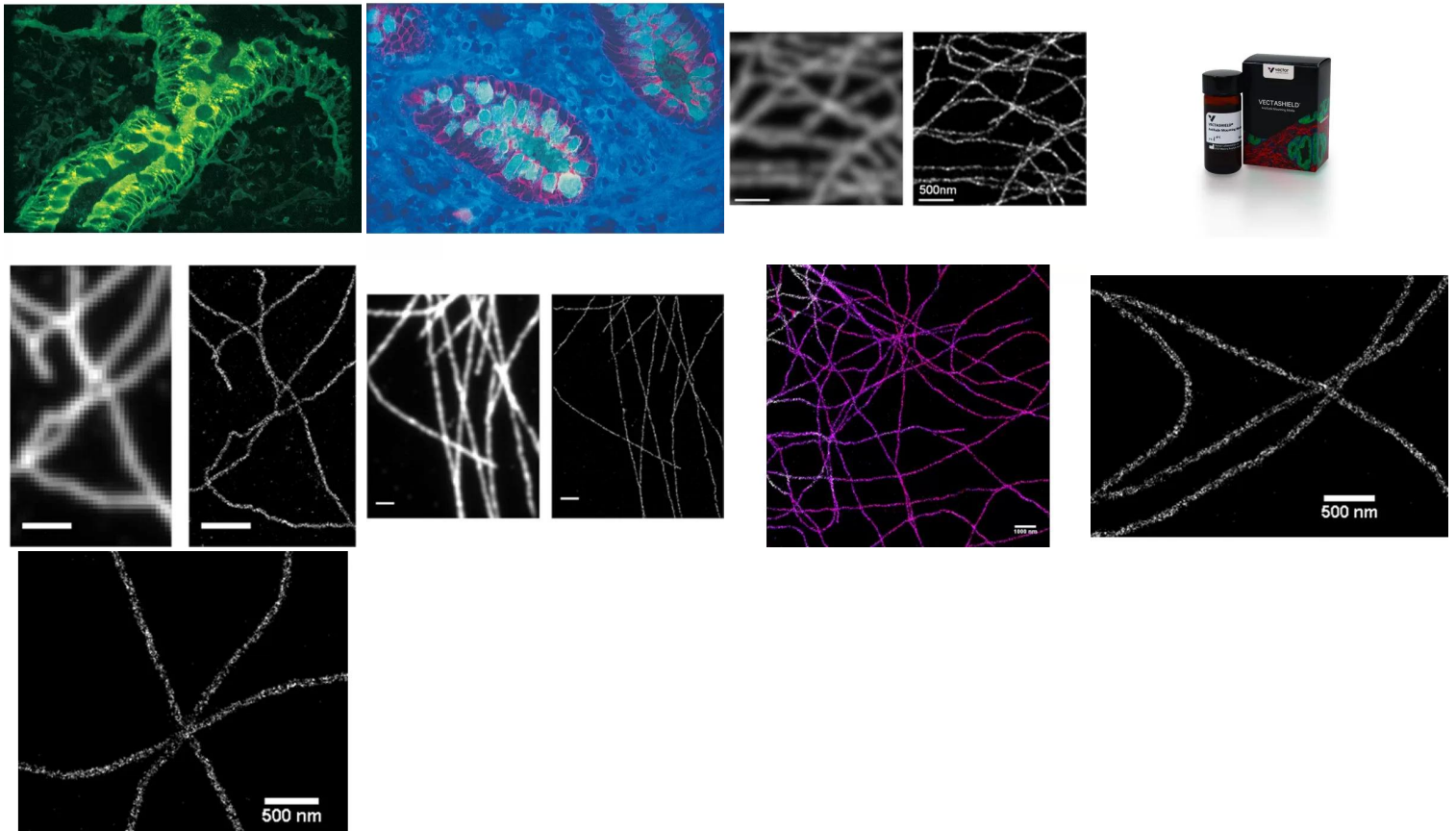

**For research use only. Not intended for therapeutic or diagnostic use in animals or humans.**